

FitRx

max RECOVER massager



USER MANUAL




INTRODUCTION & FEATURES

Relieve back, shoulder pain and tensed muscles with the FitRx massager. This ultra-quiet and portable massager is your answer to soothed and relaxed muscles. The FitRx massager is optimized to offer a personalized massage experience. Choose from six distinct interchangeable massage heads, select the speed of your massage using the adjustable frequency on the screen indicator and start the massage. Compact and portable the FitRx massager combines several massage functions in a smart way for rapid relief from fatigue, soreness and stress.

PACKAGE CONTENTS:

- FitRx Massage Gun
- Charger
- Six Massage Heads
- Manual
- Case

PRODUCT OVERVIEW & FUNCTIONS

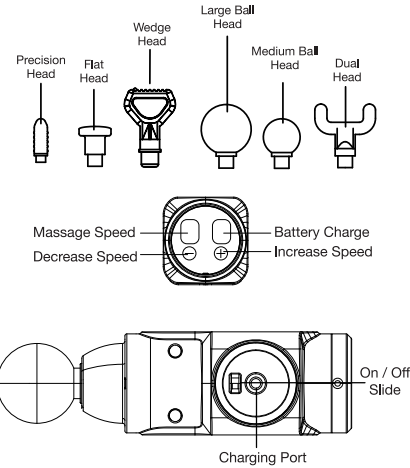

USER INSTRUCTIONS

- To turn the unit on slide the ON/OFF button at the bottom of the massager to the ON position.
- The numbers on the right side of the screen represent the battery power: 99/80/60/40/20 from full to low. The numbers on the left side 01-06 reflect the speed of the massager.
- Select the massage head you would like to use and insert it into the round hole in front of the massager.
- Press the (+) button on the bottom of the screen to start the massager
- Each press on the (+) button increases the speed of your massage.
- For the 01 speed the pulse frequency is 1200 rpm. Every additional press on the (+) button increases the speed by 400 rpm. The highest is 06 reaching a top speed of 3200 rpm.
- Press the (-) button to decrease the speed.
- A long press on (+) or (-) buttons for five seconds will automatically stop the massager.
- After completing the massage switch the battery (OFF) at the bottom of the massager and the screen will automatically turn off.

NOTE: The FitRx massager is equipped with a safety function it will stop running automatically if it has been running continuously for more than 10 minutes or if the application force is strong. Simply turn it on again if you want to keep using it.

CHARGING INSTRUCTIONS

When the battery is exhausted it needs to be fully charged for at least 1-2 hours. The battery working time is determined by the speed and frequency of the massage. While charging the power adapter LED light will glow red once fully charged it will turn green.

SPECIFICATIONS

- Working Voltage: 12V
- Battery Voltage: 12V
- Output Power: 12 Watts
- Charging Time: 1-2 Hours
- Six interchangeable massage heads
 - Large ball head / Medium ball head
 - Precision head
 - Dual head
 - Wedge head
 - Flat head

SAFEGUARDS:

To reduce the risk of electric shock, fire and injury follow these instructions carefully:

- This device is recommended for use by adults only.
- Do not use the massager directly on the surface of your skin. To avoid injury the massaged area should be covered by clean clothes or a dry towel.

- Gently press and move the massager; do not massage the same part of the body for more than 60 seconds.
- Use high speed on soft body tissue only. Do not use high speed and pressure on the head, bones and joints.
- Always choose the right massage head for the respective body part
- Frequent use of high pressure may cause abrasions. If you experience pain or discomfort stop using the massager immediately
- Keep fingers and bodily hair away from the telescopic rod of the device and the heat dissipation hole to prevent pinching and injury.
- Do not place anything in the heat sink of the motor.
- The massager should not be submerged in water. Do not allow water/liquids to enter the massager from the motor cooling hole.
- Use the massager with its original charger only.
- Do not modify or make changes to the device's equipment.
- Do not run and charge the device in an unmanned state.

DISCLAIMER: Read the instructions carefully before using the massager. The use of this equipment is subject to following its instructions closely. Tzumi is not liable for any damage or loss caused by inappropriate use of this device.

Do not use this massager without the approval of your doctor if you are: pregnant, diabetic, have a pacemaker, recovering from surgery, suffer from migraines, neuropathy disorders, disc protrusion or spinal inflammation. If you feel pain or discomfort while using the massager stop using it immediately and consult your doctor or physiotherapist.

tzumi® Distributed by tzumi inc. NY, NY, 10016. ©Copyright 2020 tzumi. All Rights Reserved. Made In China, Electronics Co., Ltd. All pictures shown are for illustration purpose only. Actual product may vary due to product enhancement. Other trademarks and trade names on this package are those of their respective owners.